




Simply *Guided*

Neoss® Guide Kit – Hard Stop



Intelligent Simplicity

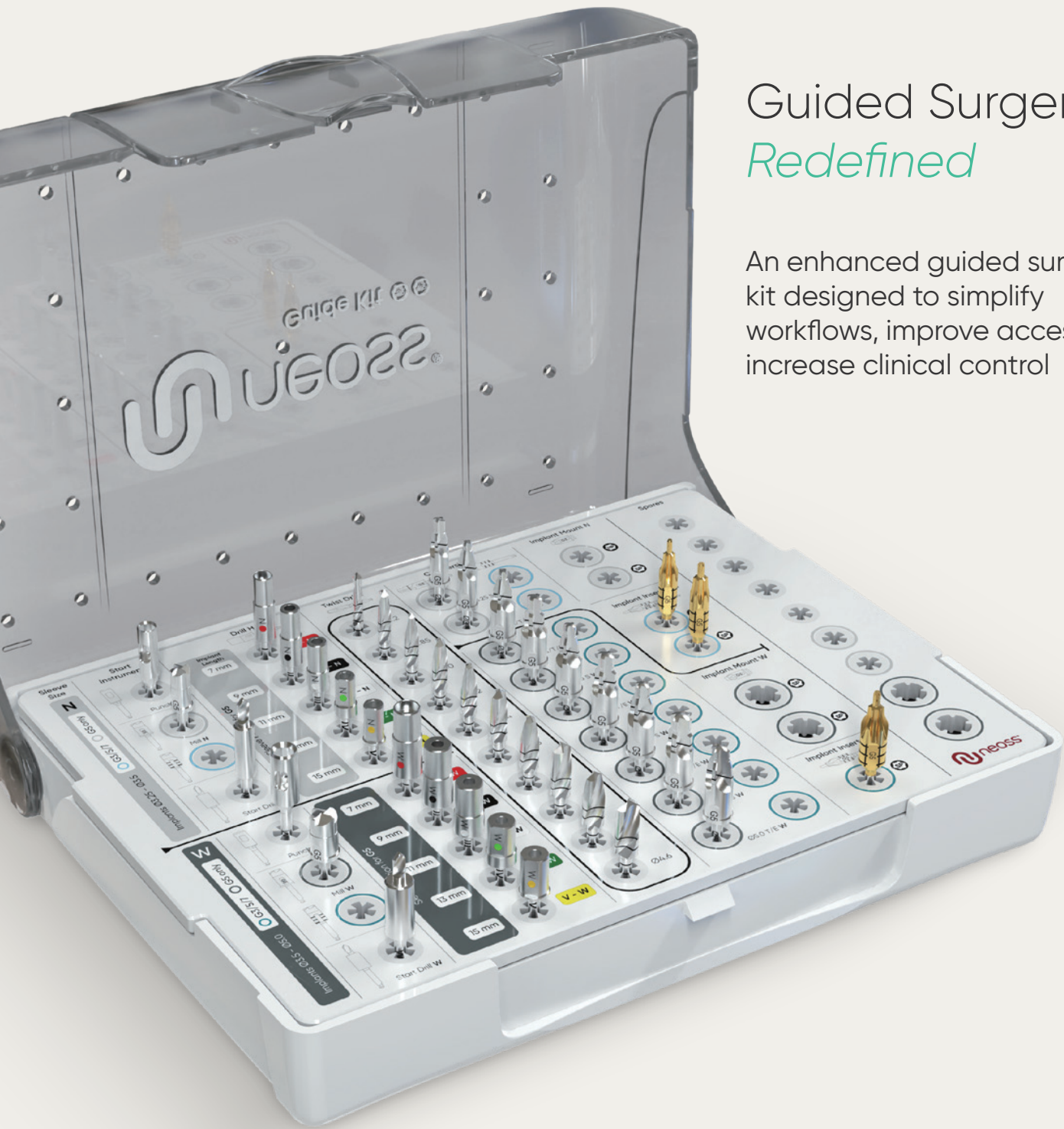

Guided Surgery, *Redefined*

An enhanced guided surgery kit designed to simplify workflows, improve access and increase clinical control

Shorter tools, *eliminate guide keys*

New hub & drill system

The drill hub eliminates the need for separate guide keys to simplify handling during drilling protocols. The overall vertical height of tools is significantly reduced for better access, particularly in the posterior region.





One guided surgery kit for all three Neoss systems
Newly designed drilling protocol

All our implant systems can be inserted with a single set of drills and a single guide kit, reducing the complexity of surgical preparation and simplifying your inventory.

Intuitive design for fast adaptation

Streamlined workflow

The Hard Stop Kit setup is designed to work with the fixed G5 (5mm) sleeve offset, with one drill hub corresponding to one implant length, and hard stops on every bone cutting instrument. This intuitive design ensures effortless depth control, reduced learning curve and improved procedural predictability.

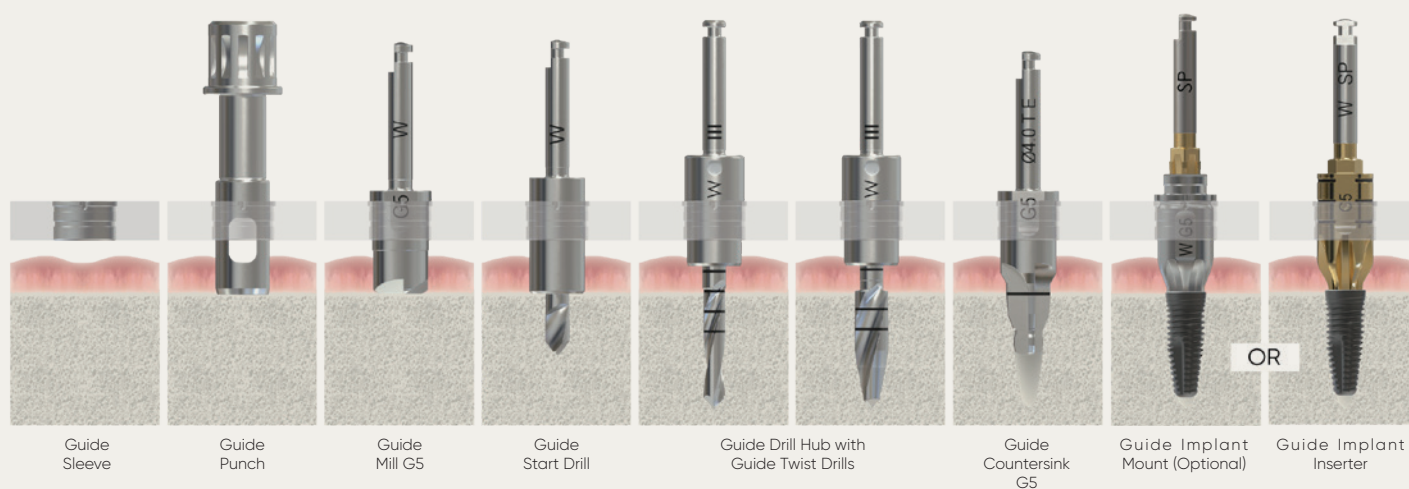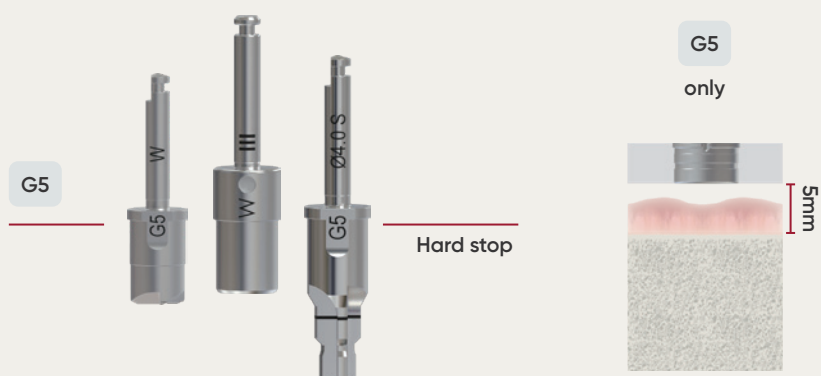

The new sleeve design enhances seating stability and supports precise orientation control. Open sleeves reduce the vertical height requirement, improving access in challenging cases.


	Closed Sleeve	Open Sleeve	Implant Diameter
Narrow – N (Ø4.2)			Ø3.25 – Ø3.5
Wide – W (Ø5.5)			Ø3.5- Ø5.0

Achieve precise *pre-designed* implant orientation

Full digital workflow compatibility

Integrated orientation control on the implant inserter and sleeves, which is also fully reflected in our digital design databases, enables accurate compatibility with pre-designed and prefabricated temporary restorations for immediate loading cases.

Design databases are available in leading planning software programs, with more on the way.

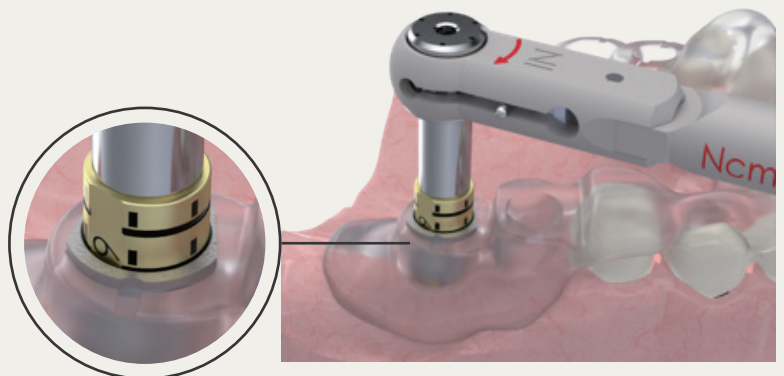

Assortment

Neoss® Guide Kit Ø3.25 – Ø5.0 Hard Stop

Art. No. Hard stop and flexible kit and ad on

51200	Neoss® Guide Kit Ø3.25 – Ø5.0 G5 Hard Stop
51260	Neoss® Guide Instrument Kit <i>(refer to next page)</i>

51200 includes:

Art. No. Description

51201	Neoss® Guide – Punch – Narrow
51246	Neoss® Guide – Mill – Narrow G5
51203	Neoss® Guide – Start Drill – Narrow
51204	Neoss® Guide – Punch – Wide
51247	Neoss® Guide – Mill – Wide G5
51206	Neoss® Guide – Start Drill – Wide

51207	Neoss® Guide – Drill Hub I, Narrow
51208	Neoss® Guide – Drill Hub II, Narrow
51209	Neoss® Guide – Drill Hub III, Narrow
51210	Neoss® Guide – Drill Hub IV, Narrow
51211	Neoss® Guide – Drill Hub V, Narrow


51212	Neoss® Guide – Drill Hub I, Wide
51213	Neoss® Guide – Drill Hub II, Wide
51214	Neoss® Guide – Drill Hub III, Wide
51215	Neoss® Guide – Drill Hub IV, Wide
51216	Neoss® Guide – Drill Hub V, Wide

51217	Neoss® Guide – Twist Drill, Ø2.2 × 7-15 mm
51218	Neoss® Guide – Twist Drill, Ø2.85 × 7-15 mm
51219	Neoss® Guide – Twist Drill, Ø3.0 × 7-15 mm
51220	Neoss® Guide – Twist Drill, Ø3.2 × 7-15 mm
51221	Neoss® Guide – Twist Drill, Ø3.4 × 7-15 mm
51222	Neoss® Guide – Twist Drill, Ø3.6 × 7-15 mm
51223	Neoss® Guide – Twist Drill, Ø3.9 × 7-15 mm
51224	Neoss® Guide – Twist Drill, Ø4.1 × 7-15 mm
51225	Neoss® Guide – Twist Drill, Ø4.4 × 7-15 mm
51226	Neoss® Guide – Twist Drill, Ø4.6 × 7-15 mm

51249	Neoss® Guide – Countersink, Ø3.25 S – Narrow G5
51250	Neoss® Guide – Countersink, Ø3.5 S/T/E – Narrow G5
51251	Neoss® Guide – Countersink, Ø3.5 S/T/E – Wide G5
51252	Neoss® Guide – Countersink, Ø4.0 S – Wide G5
51253	Neoss® Guide – Countersink, Ø4.5 S – Wide G5
51254	Neoss® Guide – Countersink, Ø5.0 S – Wide G5
51255	Neoss® Guide – Countersink, Ø4.0 T/E – Wide G5
51256	Neoss® Guide – Countersink, Ø4.5 T/E – Wide G5
51257	Neoss® Guide – Countersink, Ø5.0 T/E – Wide G5


51238	Neoss® Guide – Implant Inserter, SP, Wide
51239	Neoss® Guide – Implant Inserter, SP, Narrow
51240	Neoss® Guide – Implant Inserter, NP, Narrow

51241	Neoss® Guide Kit System Box
-------	-----------------------------

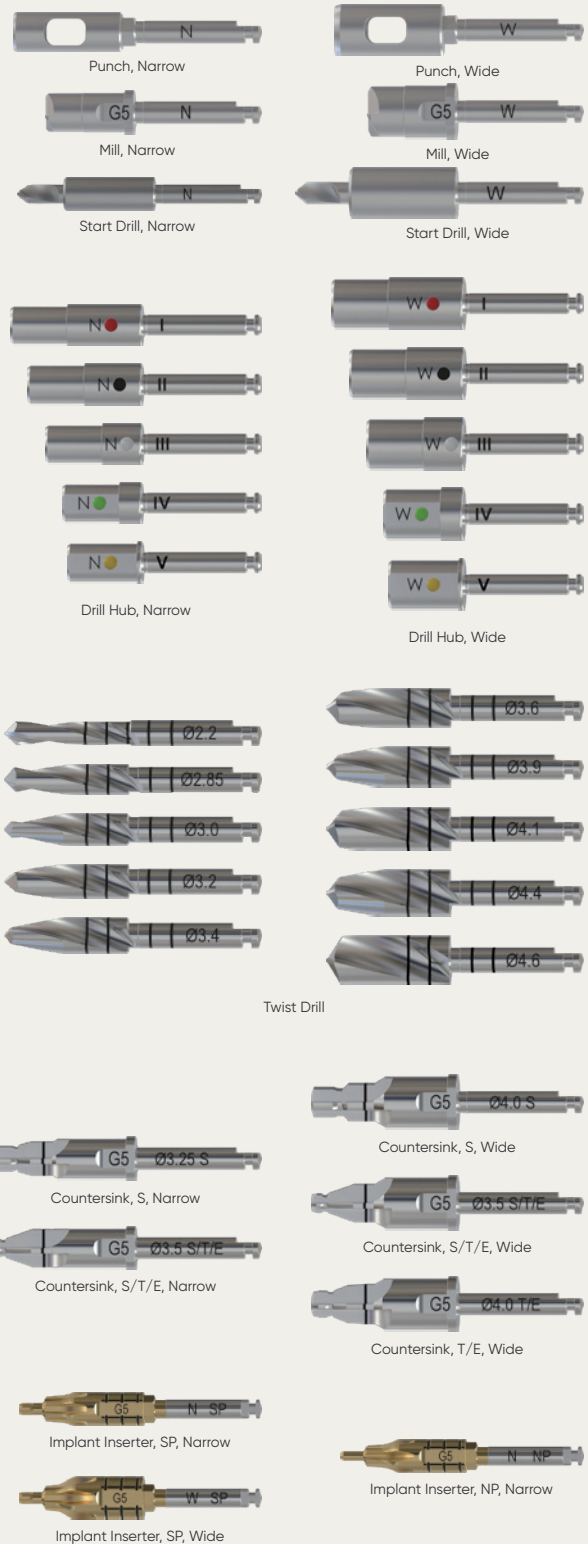

51260 includes:

Art. No.	Description
51121	Ratchet
51144	Wrench Adapter
51137	Implant inserter, SP 17mm
51126	Manual Handle
51139	Neo Screwdriver, Machine 22mm

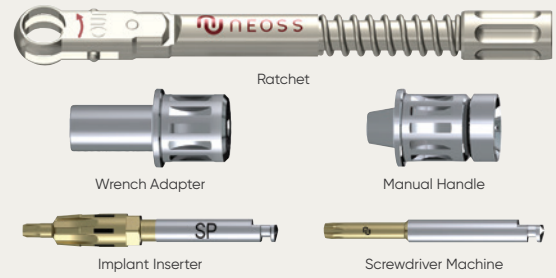

Guide Sleeves

Guide Sleeves can be ordered separately for Guide.

Art. No.	Description
51242	Neoss® Guide - Open Sleeve, Narrow - 5 pcs
51243	Neoss® Guide - Closed Sleeve, Narrow - 5 pcs
51244	Neoss® Guide - Open Sleeve, Wide - 5 pcs
51245	Neoss® Guide - Closed Sleeve, Wide - 5 pcs



Implant Mounts (Optional)

Implant Mounts are optional and can be ordered separately.

Art. No.	Description
51235	Neoss® Guide - Implant Mount, NP, Narrow G5
51236	Neoss® Guide - Implant Mount, SP, Narrow G5
51237	Neoss® Guide - Implant Mount, SP, Wide G5



You can purchase 51261 Neoss® Guide add-on-kit Countersinks and Mill G3/5/7 for extra components in flexible mode. Please contact your local Neoss Representatives for the options.

AUSTRALIA & NEW ZEALAND

Neoss Australia Pty. Ltd.
 G04 Ground Floor North Tower, John Oxley Centre
 339 Coronation Drive
 Milton Queensland 4064
 Australia
 T: + +61 7 3216 0165
 E: info.au@neoss.com

Order your products
 24/7 in our Online Stores



Australia
 at shop-au.neoss.com



New Zealand
 at shop-nz.neoss.com

For our Terms & Conditions consult your local website on www.neoss.com

Contact your local representative for product availability details.
 All products are not available in every market.

17518_0 EN 26-03 © Neoss Group, 2026. Copyrights, design rights and trademarks, Neoss documents, software and designs may not be reprinted, copied or published in whole or part, without the written authorization of Neoss Group. Neoss and the N logo are registered trademarks of Neoss Group. NeoScan™, NeoPro™, and NeoConnect™ are trademarks of Neoss Group.